

S630

DENTAL UNIT

SIMPLE AND COMFORTABLE
PRACTICAL AND RELIABLE



S630

DENTAL UNIT



ENGINEERING ELECTRIC CHAIR

Comfortable chair cushion design

The seat design is based on the human body shape and provides perfect comfort and support, especially for patients who have been treated for a long time.



S630

DENTAL UNIT



- ▶ **LED Lamp**
The handle is removable for cleaning, sterilization.



- ▶ **Sensor on/off** (white light, yellow light and mixed 3 kinds of light source mode manual conversion).
Reduces eye fatigue and precise colorimetry.



- ▶ **Assistant holder combination device** (Delay off)
Aluminum handle can be removed for cleaning and disinfection.



- ▶ **Assistant holder silicone pad**
Scratch resistant, easy to clean and disinfectable.



- ▶ **Aluminum instrument plate (Interlock)**
The handle is removable for disinfection and the hanger is rotatable.



- ▶ **Stainless steel tray**
Large size, easy to clean.



- ▶ **External water storage bottle**
One-key conversion, one-key switch, convenient and fast.



- ▶ **Removable strength and weakness suction filtration device**
Easy to clean.

S630

DENTAL UNIT



- ▶ **Double joint headrest**
The fold-out design makes the patient more comfortable and at ease during treatment.



- ▶ **Silent electric chair (Security Protection)**
The lowest chair position is 380mm design, which is convenient for children and elderly patients.




- ▶ **Rotatable handrail design**
Aluminium alloy / engineering plastic face cover.



- ▶ **Emergency stop switch**
One-click stop, safe and convenient.



- ▶ **Ceramic spittoon, automatic water supply system**
Removable / rotating, rinsing mouth ration water.

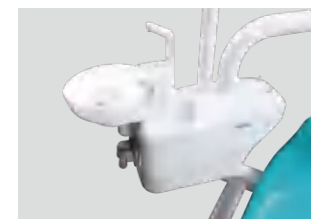

- ▶ **The chassis can be rotated 45 degree**
Easy to operate with four hands.




- ▶ **Multifunctional foot pedal**
Control up and down / supine / flushing cup / water supply / oral lights / spitting position / handpiece working water status / speed.



- ▶ **Water/air/electricity master switchgear (Anti-suction device)**
One-click on/off convenient and fast.

S630

DENTAL UNIT

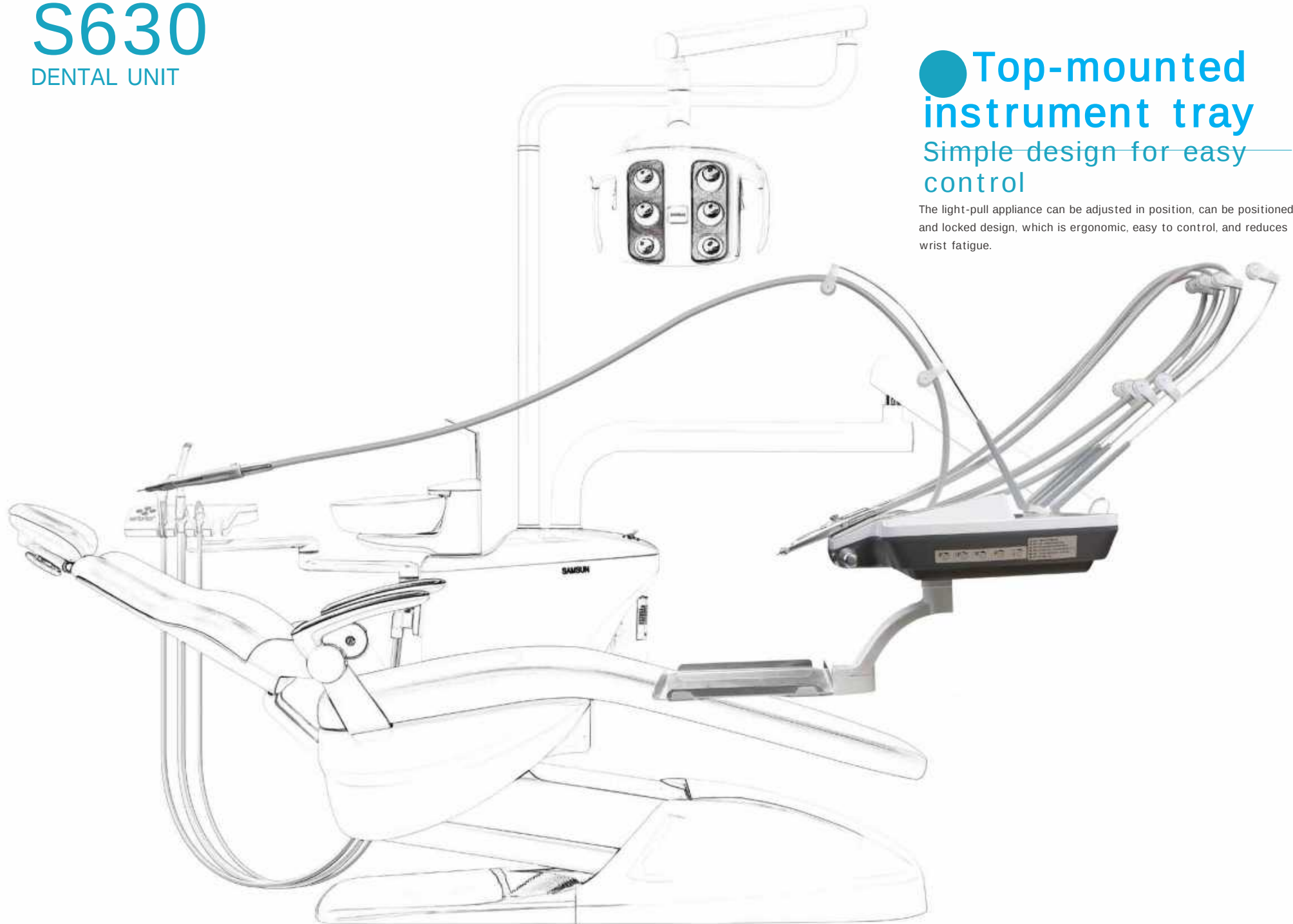

● Top-mounted instrument tray

Simple design for easy control

The light-pull appliance can be adjusted in position, can be positioned and locked design, which is ergonomic, easy to control, and reduces wrist fatigue.

S630

24V DC silent electric chair(safety protection)	■
9 programs control system(integrated circuit board)	■
Auto spittoon flushing and cup filler control system	■
Top mounted instrument tray (interlock/air break)	■
Stainless steel tray (rotatable)	■
Ceramic spittoon (rotatable/removable)	■
Assistant holder combination (folding and rotatable)	■
S2LED lamp(sensor on/off/illuminancea djustment)	■
Lamp armset(with damping)	■
Multi-function foot pedal	■
Strong/weak suction(disassembly aluminum handle)	■
Folding double joint head rest	■
Double armrests	■
Seamless seat cushion(anti-static PU leather)	■
Medical gradewater/trachea	■
X film viewer	■
USB cable(integrated line)	■
External water storage	■
Water/air/electricity main switch(anti-suction)	■
3-way syringe(hot/cold)	■
Constant temperature heating	■
Build-in floor box	■
S60 doctor stool	■

Options

S1 LED lamp	☐
S3 LED lamp	☐
Disinfectior for water	☐
Built out floor box	☐
Inductive water supply	☐
Fiber leather	☐
Glasss pittoon	☐
Intraoral camera(withbracket)	☐
Central suction machine	☐
Down mounted instrument tray	☐
Cart type instrument tray	☐
S30 doctor stool	☐
LED light curing machine (built-in)	☐
Teeth scaler(built-in)	☐
Imported/domestic high-speed handpiece	☐
Imported/domestic low-speed handpiece	☐
Imported/domestic electric motor(built-in)	☐

Optional Accessories





S1 LED LAMP S3 LED LAMP Disinfectant for water **Infrared camera**



Aluminum alloy instrument plate	Cart type instrument tray	S30 doctor stool	S60 doctor stool
---------------------------------	---------------------------	------------------	------------------

Construction reference drawings

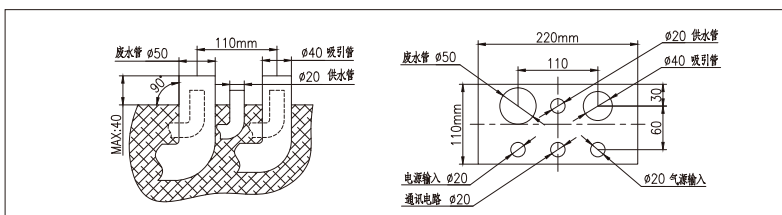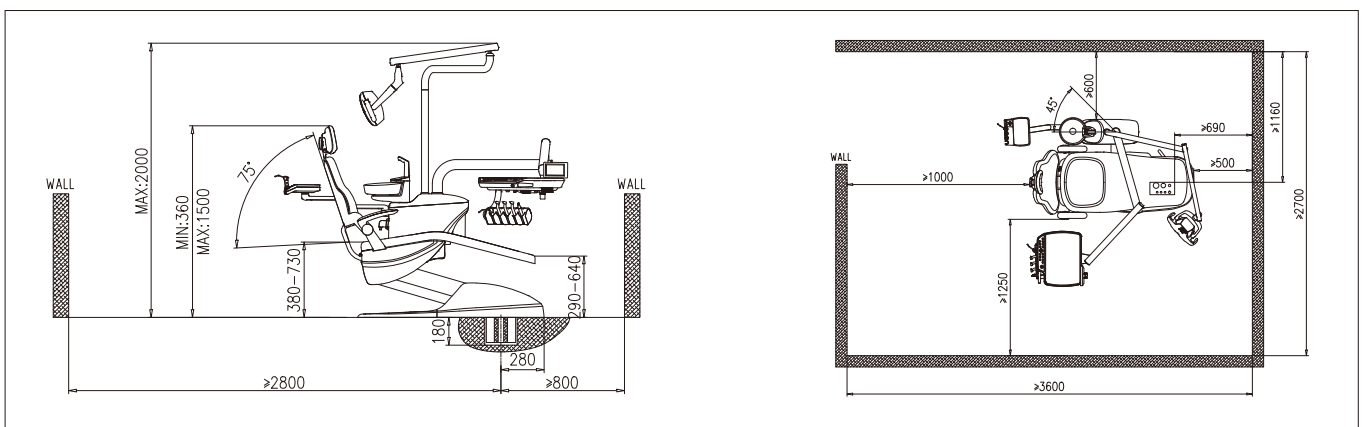

Available colours

Sewed ECO-leather cushion(optional)



H6 H7 H8 H9 H10

Seamless PU leather



Technical parameters

■ Power supply voltage	AC220V/50HZ
■ Water pressure	0.2MPa-0.4 MPa
■ Air source pressure	0.5MPa-0.8MPa
■ Air displacement	>55L/min