

S680

DENTAL UNIT

SIMPLE AND COMFORTABLE
PRACTICAL AND RELIABLE



S680

DENTAL UNIT



ENGINEERING ELECTRIC CHAIR Comfortable chair cushion design

The seat design is based on the human body shape and provides perfect comfort and support, especially for patients who have been treated for a long time.



S680

DENTAL UNIT



- ▶ **LED Lamp**
The handle is removable for cleaning, sterilization.



- ▶ **Sensor on/off** (white light, yellow light and mixed 3 kinds of light source mode manual conversion)
Reduces eye fatigue and precise colorimetry



- ▶ **Assistant holder combination (off delay)**
Aluminum handle can be removed for cleaning and disinfection.



- ▶ **Assistant holder silicone pad**
Scratch resistant, easy to clean and disinfectable




- ▶ **Aluminum instrument plate (interlock)**
The handle is removable for disinfection and the hanger is rotatable



- ▶ **Stainless steel tray**
Large size and easy to clean



- ▶ **External water storage bottle**
One-key conversion, one-key switch, convenient and fast

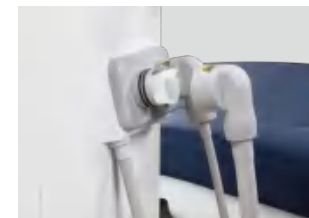

- ▶ **Removable strength and weakness suction**
filtration device easy to clean

S680

DENTAL UNIT



- ▶ **Double joint headrest**
The fold-out design makes the patient more comfortable and at ease during treatment



- ▶ **Silent electric chair (Security Protection)**
The lowest chair position is 380mm design, which is convenient for children and elderly patients



- ▶ **Rotatable handrail design**
Aluminium alloy/engineering plastic face cover



- ▶ **Emergency stop switch**
One-click stop, safe and convenient



- ▶ **Ceramic spittoon, automatic water supply system**
Removable/rotating, rinsing mouth ration water.



- ▶ **Floor-standing box**
Installed on the ground, high stability



- ▶ **Multifunctional foot pedal**
cup/water supply/oral lights/ spitting position/handpiece working water
Control up and down /supine/ flushing



- ▶ **Water/air/electricity switchgear (Anti-suction device)**
One-click on/off convenient and fast

S680

24V DC silent electric chair(safety protection)	■
9 programs control system(integrated circuit board)	■
Auto spittoon flushing and cup filler control system	■
Floor-standing box	■
Cart type instrument tray (interlock/air break)	■
Stainless steeltray(rotatable)	■
Ceramic spittoon (rotatable/removable)	■
Assistant holder combination (folding and rotatable)	■
S3 LED lamp(sensor on/off/illuminancea djustment)	■
Lamp armset(with damping)	■
Multi-function foot pedal	■
Strong/weak suction(disassembly aluminum handle)	■
Folding double joint head rest	■
Seamless seat cushion(anti-static PU leather)	■
Double armrests	■
Medical gradewater/trachea	■
X film viewer	■
USB cable(integrated line))	■
External water storage	■
Water/air/electricity main switch(anti-suction)	■
3-way syringe(hot/cold)	■
Constant temperature heating	■
Build-in floor box	■
S60 doctor stool	■
S30 doctor stool	■

Options

S1 LED LAMP	☐
S2 LED LAMP	☐
Disinfector for water	☐
Built out floorbox	☐
Inductive water supply	☐
Glassspittoon	☐
Intraoral camera(with bracket))	☐
Central suction machine	☐
Down-mounted instrument tray	☐
Cart type instrument tray	☐
Teeth scaler machine(built-in)	☐
LED light curing machine(built-in)	☐
Fiber leather	☐
Imported/domestic high-speed handpiece	☐
Imported/domestic low-speed handpiece	☐
Imported/domestic electric motor(built-in)	☐

Optional Accessories





S1 LED LAMP S2 LED LAMP Disinfector for water Intraoral camera



Top-mounted instrument tray	Aluminum alloy instrument plate	S30doctor stool	S60doctor stool
-----------------------------	---------------------------------	-----------------	-----------------

Construction reference drawings

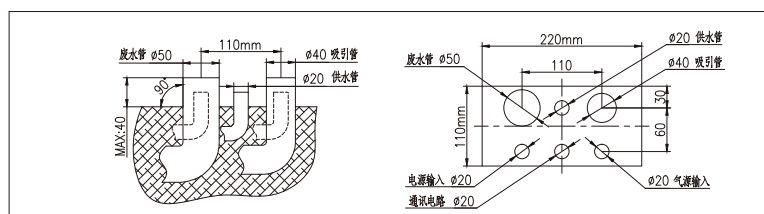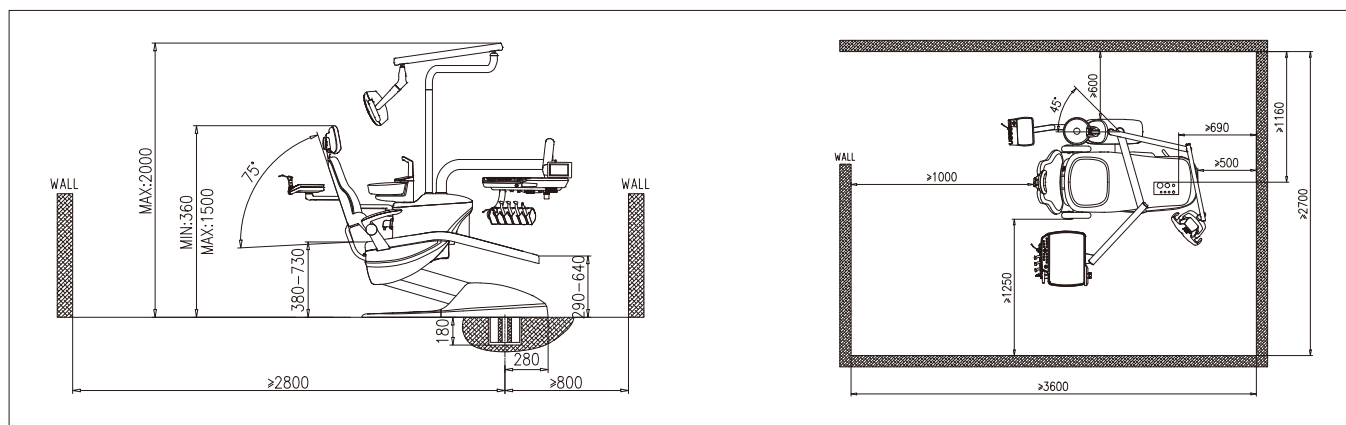


Available color

Sewed ECO-leather cushion(optional)



Seamless PU leather

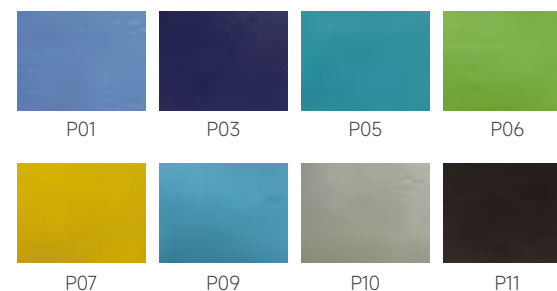

Technical parameters

■ Power supply voltage	AC220V/50HZ
■ Water pressure	0.2MPa-0.4MPa
■ Air source pressure	0.5MPa-0.8MPa
■ Air displacement	>55L/min